

Now with the updated
CryoCube® F101h



Pure. Ultracold.

The Eppendorf Ultra-Low Temperature Freezer Family



»Design the freezer around the sample.«

After nearly 40 years, our freezer journey still continues. By combining the longevity and quality of our existing freezers with future-proven sample monitoring and management systems, we have designed an icon for -86 °C ultra-low temperature freezers in the 700 L class, the green CryoCube® F740 freezer series. But that is just the beginning of the journey.



Concerned about your sample safety?

- > Temperature monitoring and data export for documentation*
- > Controlled access to samples for higher security*
- > Dedicated alarm and backup system for 24/7 sample safety

> See page 6 for more information


* CryoCube F740 » series models.



Suffering from energy bills?

- > High-efficiency insulation for low power consumption
- > Special gaskets on inner and outer doors to keep the cold inside
- > High-efficiency compressor systems provide excellent performance with minimized energy usage/carbon footprint

> See page 10 for more information

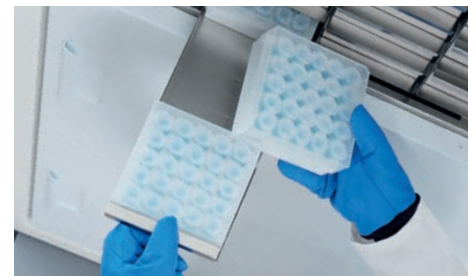

Lost in samples?

- > Broad range of metal racks for flexible vessel usage
- > Storage boxes that fit standard vessel sizes
- > Sample management system to keep your inventory in order
- > Broad range of vessel types including plates, microtubes, conical tubes, and 2D barcoded vessels


> See page 17 for more information

How Much Effort Did You Spend on Your Sample?

Did you ever calculate the value within your freezer? All the hours, days, weeks, and months spent on your samples? More than 50,000 samples fit into large, standard ULTs. The value of every single sample differs – from simple buffers to high-value cell extracts, expensive enzymes, or very rare sample material. Assuming an average value of 10 Euro or Dollar per vial, the total value already reaches 500,000 Euro or Dollar.

When predominantly storing high-value samples which are a result of months of work – you reach a value in your freezer far beyond half a million Euro or Dollar.


Your new ULT freezer is not limited to being a storage room for samples. This instrument is your assurance for long-term scientific success, it guards the results of your work.

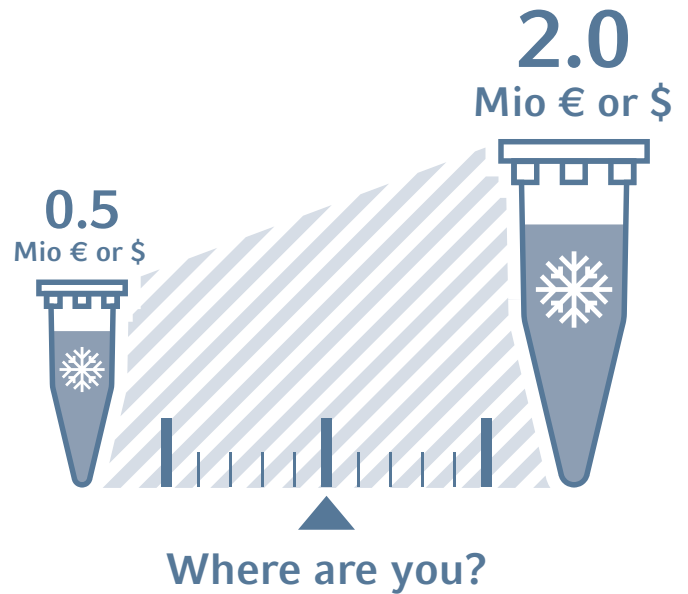

Life Science Revolves Around the Sample

Your processes and workflows require you to invest a lot of time and a lot of resources. So much of your success depends on the care and handling of your valuable samples.

After spending countless hours considering the many devices and steps utilized to create the sample, have you adequately considered where you will store it? Does your long-term storage solution match the value of your processes and samples?

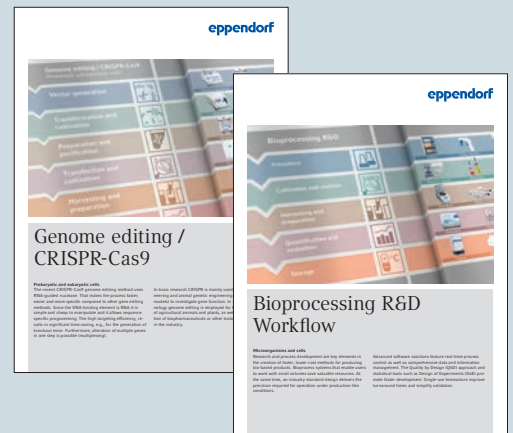

Remember the last time it happened?


»...it was 2 a.m. when the guard called me: Sir, there is an alarm in your lab, you'd better come over...«

»...it's a nightmare, you enter the lab, the water is already dribbling down the freezer door, and you know – it's all gone...«

»...last summer our freezer failed, we lost 2 years of work and 650,000 Euro...«



After all this work, do you know where your samples are?

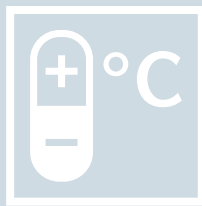

Bioprocessing Workflow

Key & benefits:
The standard method for cell culture and bioprocessing is the use of flasks and bioreactors. However, the use of Eppendorf's Bioprocessing Workflow offers a more efficient and cost-effective solution. It allows for the production of high-quality cell cultures and bioprocesses in a controlled and reproducible manner. The workflow is designed to be flexible and scalable, allowing for the production of small-scale cultures for research and development, as well as large-scale cultures for commercial production.

Stem Cells Workflow

Research and applications:
The stem cell workflow is a key technology in regenerative medicine and cell-based therapies. It allows for the generation of pluripotent stem cells (PSCs) from somatic cells, which can then be differentiated into various cell types for research and clinical applications. The workflow is designed to be efficient and reproducible, allowing for the production of high-quality stem cells for research and clinical use.

Next Generation Sequencing

Eppendorf NGS workflow efficiency:
The workflow is designed to be efficient and reproducible, allowing for the production of high-quality NGS libraries. The workflow is designed to be flexible and scalable, allowing for the production of small-scale libraries for research and development, as well as large-scale libraries for commercial production.

Molecular Biology Workflow

Key & benefits:
The molecular biology workflow is a key technology in the study of gene expression and regulation. It allows for the production of high-quality RNA libraries for sequencing and analysis. The workflow is designed to be efficient and reproducible, allowing for the production of high-quality RNA libraries for research and clinical use.

Cell Biology Workflow

Key & benefits:
The cell biology workflow is a key technology in the study of cell growth and differentiation. It allows for the production of high-quality cell cultures for research and clinical use. The workflow is designed to be efficient and reproducible, allowing for the production of high-quality cell cultures for research and clinical use.



Complete Control


Worried about the safety of your high value samples within your ULT freezer? With the Eppendorf ULT monitoring concept available on CryoCube F740hi and CryoCube F740hiw models, you can check the freezer performance whenever you want. All temperature data and all events are stored in the freezer control unit.

Focus on the Essentials*


The handling experience is based on our PhysioCare Concept® combined with extensive user test sessions.

- > On-board data storage means you are not dependent on auxiliary and secondary devices to keep track of freezer data
- > Easily transfer your information to your computer or colleague using the export feature and front-mounted USB port
- > Quickly and easily adjust your setpoints with a single fingertip, even while wearing gloves
- > Conveniently see all your important settings and data on a single screen with customizable views
- > Event and alarm tracking for controlled environments
- > Electronic lock system tied to user management for controlled access and higher sample safety
- > Adaptable alarm settings for individual needs
- > VisioNize® touch enabled devices: "i" series models

* CryoCube F740 "i" series models.

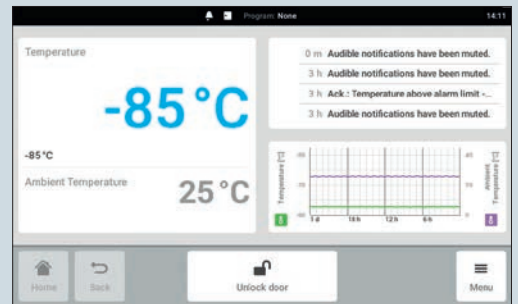

Quickly review your important parameters on the home screen*

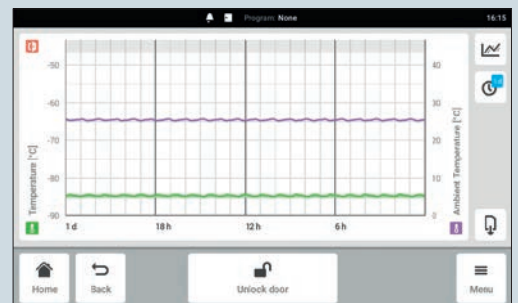

Easily review your freezer history on the temperature monitoring screen*




Proper Placement

The built-in temperature probe detects the actual temperature within the freezer. Comparing this measured value with the set value, the instrument decides if an up or down adjustment is needed.

The positioning of this probe is crucial for accurate freezer uniformity and optimized operation. Improper placement by the manufacturer can lead to inefficient cycling and a gap between your setpoint and the actual temperature inside the freezer.

Eppendorf positions the probe in the middle of the ULT chamber to receive an average value of the temperature distribution. This helps prevent inaccuracy and inefficiency; it grants peace of mind knowing your setpoint is being maintained accurately.

External temperature validation can be performed with VisioNize sense. Learn more about the VisioNize range of solutions at www.eppendorf.com/VisioNize.

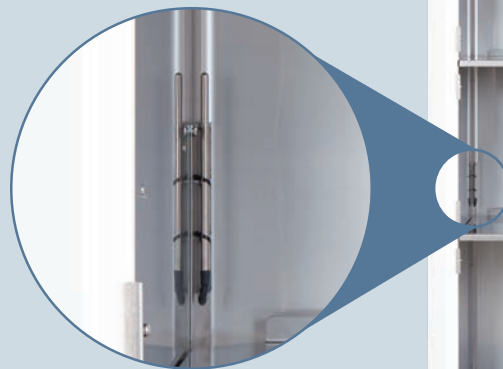

High-quality temperature probe placed in the middle of the chamber for increased accuracy and sample safety



Trust in Confirmation

Technical parameters of new ULT freezers need to be validated. On-going checks during the development phase are performed inhouse, by standard. For the final test, we at Eppendorf believe in an independent, external counter-control.

- > External test house with highly trained engineers
- > 3 units of serial production level are measured
- > Testing with defined and standardized conditions, e.g. specific room temperature, controlled humidity, validated & certificated testing instruments, and exact positioning of temperature probes
- > Documentation of all conditions and results
- > Test data are used for Eppendorf-published Performance Data Sheets





Confidence in Quality

Eppendorf has a long-standing history of innovation and quality, and our freezer manufacturing process is no exception. Each Eppendorf ULT freezer is thoroughly inspected to meet our rigorous quality guidelines. The final inspection takes more than 24 h and is based on manual and visual checks. This process is documented by an individual Certificate of Quality, complete with serial number, provided as standard for your documentation.

As your expert partner for reliable, high-quality, and efficient storage solutions, Eppendorf finely engineers each freezer with care, using only the highest quality materials to ensure a long lasting freezer that offers peace of mind and security. Eppendorf freezers are designed to offer a longterm value over the lifetime of its ownership. You can entrust your valuable samples to our ULT freezers, and you can trust your investment in all finely-engineered Eppendorf instruments.

The Future of Smart Storage



Remember Your Last 2 a.m. Wake-up Call?

Smart Storage with CryoCube & VisioNize® Lab Suite

»There is an alarm in your lab, you'd better come over...«

Due to the 50,000 valuable samples in your freezer, you rushed to the lab. Finally, the temperature was still ok, false alarm, but you were done for the day.

Just imagine, you get an alarm notification from VisioNize by email about your freezer, and in the next step, you log into the system and check the temperature graph and the events of your freezer online – all from home. Then, you can decide if you have to go to the lab immediately, if you still have some hours, or if it is another false alarm.

Your valuable samples will be now on »cloud 9« due to easy access to temperature data and events on your freezer. Your freezer becomes a Smart Storage Solution with VisioNize keeping your valuable samples safe. VisioNize offers you the possibility to remotely monitor your freezer and gives you the opportunity to set customizable notifications to meet your needs.

More peace of mind for You and Your Lab with VisioNize Lab Suite.

Monitoring



Notifications



> www.eppendorf.com/smartstorage


What Do You Need for Safe Sample Storage?

Pull-down time to -80 °C

After installation, as well as after de-frosting of the ULT freezer, the freezer needs to return to -80 °C as soon as possible. In general, a good ULT freezer needs just 4 to 5 hours to pull-down from room temperature to -80 °C, saving precious time and protecting your samples. Backup freezers are required to be ready in the shortest time possible. Based on the fast pull-down times of Eppendorf ULT freezers, backup freezers do not need to run all the time. They are ready within a short time frame. In addition, the samples within the broken unit are safe as the insulation quality supports a very slow warm-up.

Warm-up time

In cases of power failure or device break-down, the freezer insulation and the door seals are your last lines of defense against losing your valuable samples. Eppendorf ULT freezers can maintain temperatures of -50 °C or colder within your freezer for at least 8 – 10 hours*. Enough time to organize a new safe location for your samples. You do not need to hit your lab at 3 a.m..

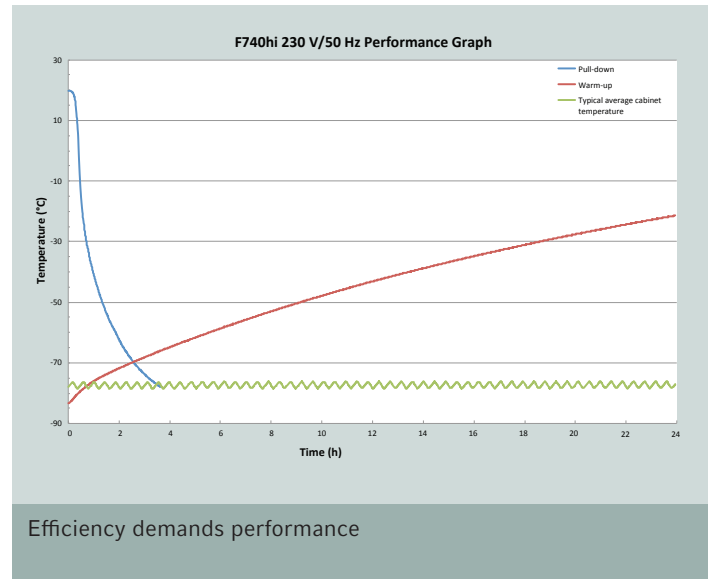

*2/3 filled

Sample integrity after door opening

How often do you open the ULT freezer door per day? For how long? It is best practice to always open and close the freezer as quickly as possible when storing new samples or attempting to locate a sample which is stored out of sight. However, this takes time.

The more time it takes, the more the temperature of both the cabinet and your frozen samples increases. The longer the door remains open, the longer it takes for the freezer to regain the set temperature. 30 to 60 seconds is a realistic time frame to add a new sample or remove a stored one.

The Door open recovery time (DOR) is the time needed to be back at -80 °C.

Door Open Recovery Time (Freezer set to -80 °C) to be back at -80 °C

	Door opening 15 seconds	Door opening 30 seconds	Door opening 60 seconds
CryoCube F570h, 230 V/ 50 Hz	14 min	19 min	27 min
CryoCube F740hi, 230 V/ 50 Hz	18 min	25 min	27 min

Safety that exceeds expectations

Besides freezer monitoring and high performance, some background safety is mandatory:


- > Secured power plug for performance safety
- > Key-lock for power switch & battery switch to avoid accidental power-off
- > Battery for alarms after power failure
- > CryoCube F740 series has a built-in voltage inverter for 100 V, 115 V, 208 V, and 230 V for high flexibility
- > Up to three different interfaces (Ethernet, RS485, BMS) for flexible connections



Cool, Cooler, Green

The coolant or cooling liquid within a ULT freezer is used to reduce and to regulate the temperature of the chamber. To reduce global warming, there is a clear tendency towards reducing classic cooling liquids as R508B and R404A.

The green or natural gases (e.g. R170, R290) are future-proof and have a very low Global Warming Potential (GWP). This is also supported by European Union Regulation (EU_517/2014).

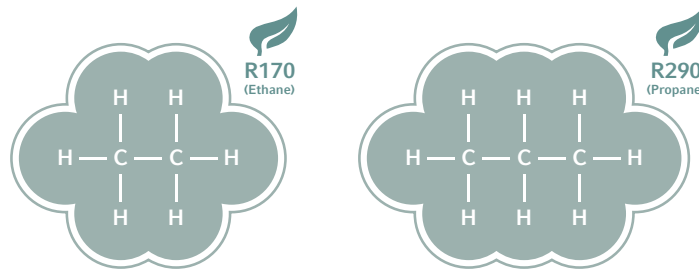

Green cooling liquids = HC gases = Natural gases = Future-Proof



Hydrocarbon cooling by Eppendorf

Remember 2008?

As one of the very first adopters of green gases at -86 °C, Eppendorf now has more than 10 years of experience in R&D, production, logistics, and service in the field.


Safety assessment for green gases:

Based on IEC 60335-2-89, there are no additional required safety instructions for using natural gases in ULT freezers as long as the cooling liquid loops are hermetically sealed and the amount per cooling liquid is limited to 150 g. All Eppendorf ULT freezers using green gases fulfill these conditions for safe usage.

The core of cooling – the compressors

Like all major suppliers of ULT freezers, Eppendorf trusts in 2-stage compressor systems. Within this 2-stage or 2-loop systems, the system is cooled down to -45 °C in the first stage and then, during the second stage, cooled down from -45 °C to -86 °C. This enables smooth and reliable operation of the individual compressors as they work in a comfort zone.

When using just one compressor, or even a Stirling engine, to cool down from +20 °C room temperature to -86 °C, more than 100 °C must be spanned in a single step. The 2-stage system is a more durable approach, supported in Eppendorf ULTs by high-value compressors from a benchmark supplier.

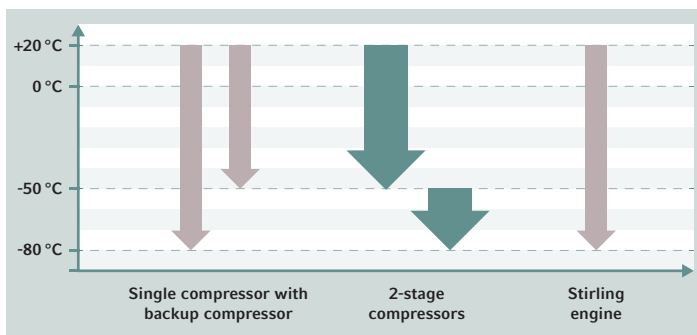


Paying Your Own Power Bill?

Even environmentally-friendly and energy-efficient ULT freezers still consume a significant amount of energy as they maintain extremely low temperatures 24/7. With today's high energy costs and focus on the environment, energy conservation has become even more important in the lab. Eco-friendly Eppendorf ULT freezers are designed to help you save energy and reduce your carbon footprint without compromising sample security.

Factors influencing the energy consumption

- > Thickness of insulation
- > Type of insulation
- > Type of cooling liquid
- > Efficiency of compressors
- > Control of compressors
- > Pressure in the cooling system
- > Style of cooling loops (diameter, length, density, etc.)
- > Seals of outer door
- > Seals of inner doors
- > Insulation of outer door
- > Insulation of inner doors

Eppendorf ULT freezer wall concepts

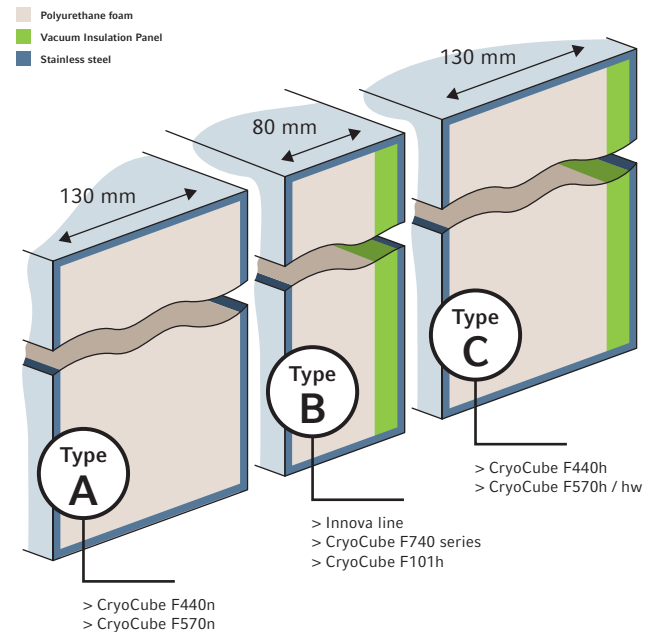

Eppendorf offers three varieties of insulation design. You can choose the design that meets your energy-saving, capacity, and cost parameters.

Keep the Cold Inside



Protect your insulation:
Displays are attached on top of insulated outer door – no heat input & no reduction of insulation




Keep the cold inside:
Efficient flat and flexible seal at outer door

Easy maintenance:
Flat seals can be quickly and easily wiped and cleaned



Double your safety:
Gaskets on insulated inner doors come standard to reduce air leakage

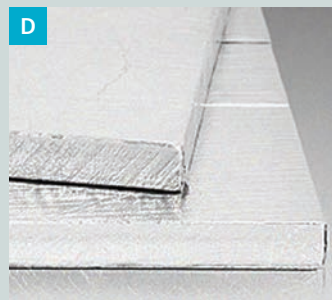


High capacity & low thermal conductivity:
Vacuum Insulation Panels provide significantly better insulation than polyurethane foam, increased capacity, as well as longterm stability

Our Part for Sustainability

- > Production moved into new facilities with state-of-the-art insulation to save on heating energy and power
 - > 100 % renewable energy in freezer facility
 - > Heat output during final individual unit testing discharges into the building heating system
 - > Plastic and cardboard compactors introduced in concert waste rebate system in place
 - > Compliance with local requirements for waste management (collection of electronic waste like circuit boards, recycling via authorized recycling organizations) and internal guidelines
 - > Usage of green cooling liquids in Eppendorf ULT freezers since 2008
 - > Green insulation foam production
 - > Packaging parts have increasing share of recycled origin
 - > Packaging is 100% recyclable
 - > ENERGY STAR® partner
 - > ACT/MyGreenLab partner
- > Marketing material is printed on certified paper with a focus on electronic mediums whenever possible
 - > Shipment by cargo ship to reduce CO₂ footprint
 - > Long lifetime of ULT freezers



Power Consumption

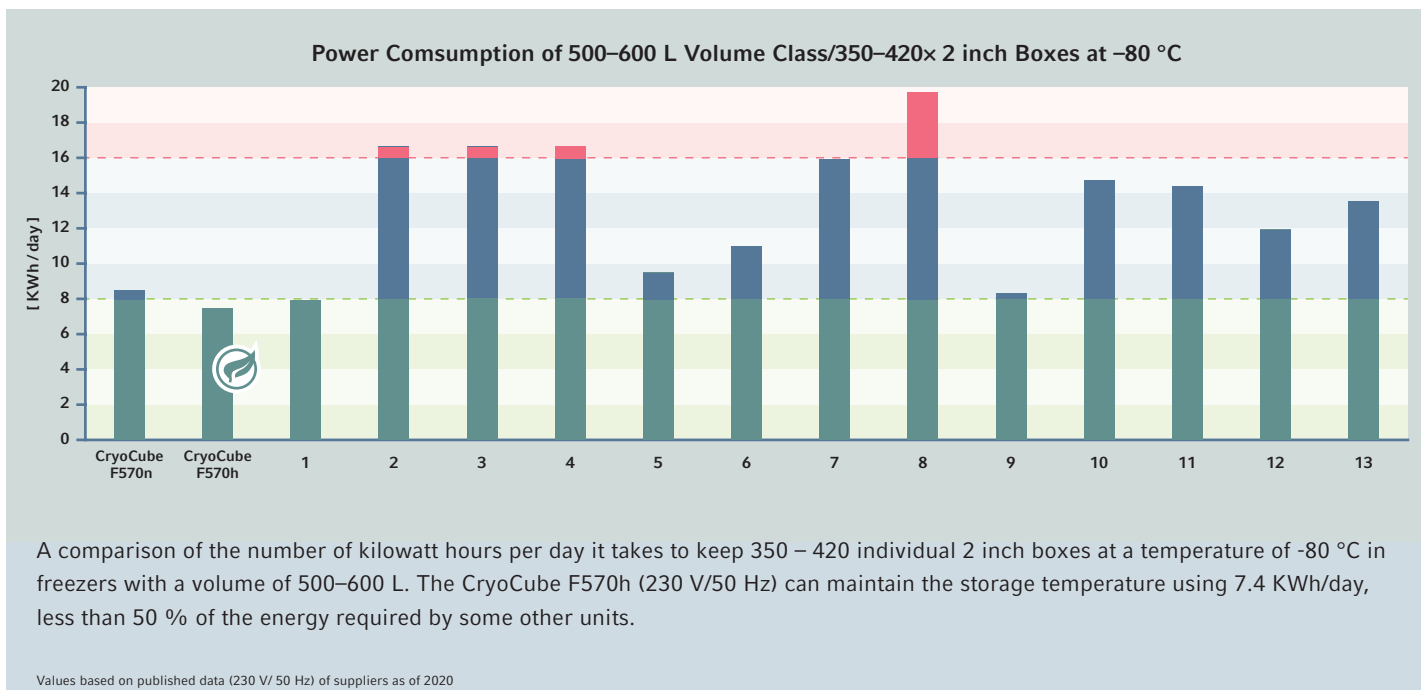

> www.eppendorf.com/Sustainability




Footprint versus efficiency

Energy-wise, chest freezers are the most economical type of freezers. Since they have sidewall insulation and are opened from the top, chest freezers can maintain their cold temperatures very well. Opening from the top reduces energy consumption per volume, especially when opening and closing the outer door of the unit on a regular basis. Less energy is lost when opening the outer door of a chest freezer than of an upright freezer as no cold can literally fall out of the freezer.

A far lower but more long-term loss of cold air is the door gasket – which is always the weakest point of a ULT freezer's insulation. As chest freezers have no vertical opening, no cold can leave the freezer via the gaskets as in an upright freezer.

However, a chest model has a much bigger footprint than an upright freezer. Transport into the building and within the building itself is more complicated for a chest freezer. Its racks (so-called "tower racks") need to be lifted up whenever you need to get a sample out of the freezer.

Air-Conditioning At the Limit? Go For Water-Cooling

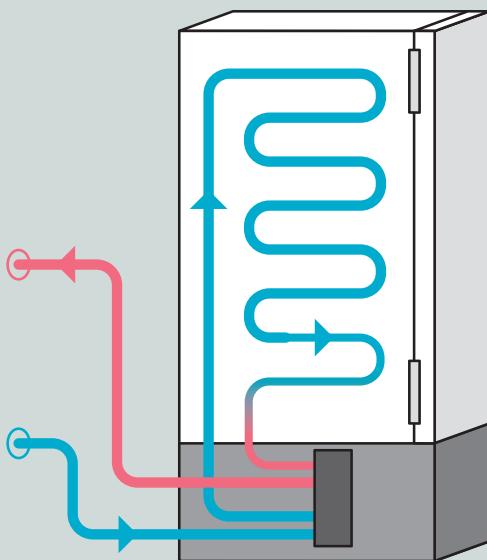

You can improve the conditions and comfort in your lab while saving energy by installing a water-cooled Eppendorf freezer. Many ULT freezers use the air in the laboratory to extract heat energy from the freezer's heat exchanger. Another option is to use a water-cooled ULT freezer connected to a facility's recirculating water system. A constant stream of water removes the heat from the heat exchanger, rather than the ambient air in the lab.


- > Reduce air conditioning power consumption by letting the water carry away the heat
- > The heated cooling water can be reused for other heat demanding systems in the facility
- > System specifications:
 - Minimum flow: 29 L/h
 - Inlet pressure: 1 – 10 bar
 - Temperature range of water: 7 °C to 25 °C

Learn more in White Paper 53:
[Air-Cooling or Water-Cooling?](#)



Eppendorf PhysioCare Concept®

The ideal Eppendorf product fulfilling the PhysioCare concept provides a comprehensive approach for the user. Starting with an ergonomic design of the product itself (e.g. smooth shape, low weight, low needed forces, intuitive user interface, etc.), fitting into an ergonomic work space, and supporting an optimal workflow within the lab.

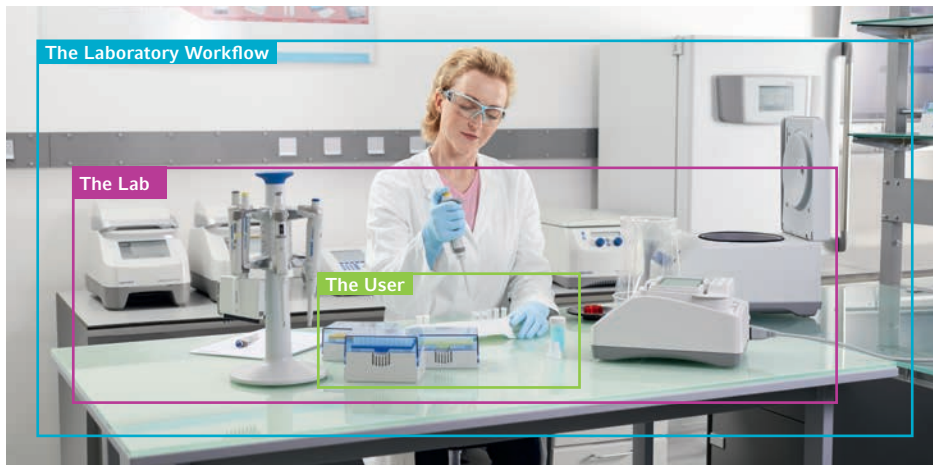

Sphere 1 – The User

Addressing the direct interaction between you and the product. Ergonomic design and product alignment optimized to the needs of the individual.

Sphere 2 - The Lab

Addressing the interaction at your workspace, resulting in ease and comfort within your workflow.

Sphere 3 - The Laboratory Workflow

General support to enhance processes around the lab and the whole company or institute.



Ergonomics? For a freezer?

Easy-to-clean:
Stainless steel interior



Enjoy opening:
Redesigned handle for easy handling, with a lock for access control and a padlock option^{*2}




Smart collection:
All cable interfaces located at one corner and 3x alarm connections for convenient access^{*1}



Be comfortable:
Eye-level automated vent-port for easy re-access with eye-level display & USB-port^{*1}



Easy access to filter:
Just flip down the cover to clean or replace^{*1}

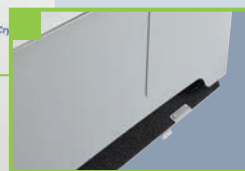

Convenient working environment:
Low noise level^{*1}



Model	Noise level (dB)A
F740h	47.8
F740 "i" series models	41.3


^{*1} CryoCube F740hi/hiw

^{*2} CryoCube F740/F570/F440 series models



In Need of Housekeeping?

The Eppendorf Storage Boxes are a complete system solution for sample storage. The comprehensive modular concept allows you to select the optimal box combinations for individual storage demands. The outer footprint dimensions (133 mm) are compatible with common freezer rack systems and the different inner grid variants provide a perfect fit for all typical tube formats: cryogenic tubes, microcentrifuge tubes, conical tubes 15, 25 and 50 mL, and other laboratory vessels. Optimize your storage and archiving of samples.

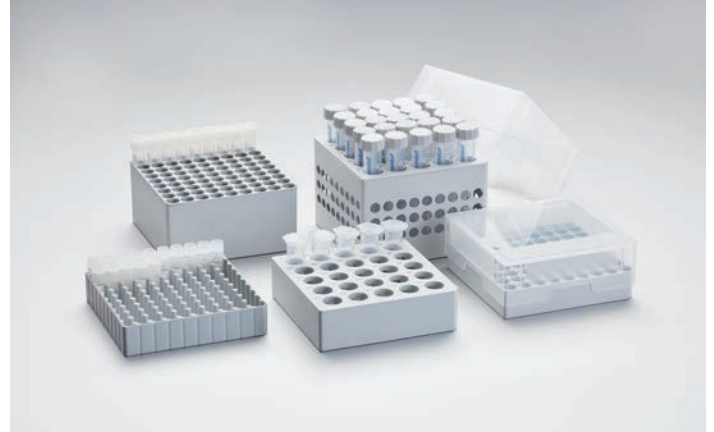

Transparent lid
for easy and fast sample inspection

Autoclavable (121 °C, 20 min)
for sterilization

5 additional inscription labels on the lid enable user-specific labeling

15 mL and 50 mL conical tube boxes with ventilation slots for an evenly distributed freezing process. Prevent tube damage or breakage

Optimal use of freezer space through flexible combination of the different formats

Made of polypropylene (PP) for high stability in freezing applications and a smooth opening and closing

Laser-labeled, high-contrast, permanent alphanumeric markings of each location for easy sample reference and minimized risk of sample mix-up

Dishwasher safe for easy cleaning

> Additional information, details on dimensions, and order numbers available at:
www.eppendorf.com/storageboxes

Burdened by Unsorted Boxes Cluttering Your Freezer?

Make your scientific life easier: store your freezer storage boxes within organized freezer racks. These metal racks optimize your freezer space and make it easy to find the items you desire. Don't waste time sifting through unsorted samples.

Drawer or side access

There are two rack styles available to make your freezer usage the most efficient. When using racks with drawers, only the relevant drawer with the freezer boxes of interest is pulled-out. All other boxes safely stay in the cold interior of the ULT freezer. Eppendorf drawers are engineered to deliver smooth movement at -80 °C for ergonomic handling. When using side access racks, accessing samples means removing a good portion of the storage boxes from the cold interior of the ULT freezer. The side access rack must be completely removed for sample access. Racks for chest freezers can be equipped with locking rods to secure freezer boxes when lifting the rack tower.

Sample storage rack selection

Though you often only hear about the safety of your samples, the method of sample storage and organization can be just as important. To select the right set of racks for optimal usage of your freezer, we advise you to consider a few questions:

1. What type of freezer do you need to equip?
2. Do you prefer stainless steel or aluminum racks?
3. Do you need drawer or side access?
4. What type of vessels do you use?
5. How can you organize your freezer in an easy and logical way?



Freedom to store: Stainless steel racks for freezer boxes of up to 136 x 136 mm footprint

Stainless steel or aluminum

By standard, there are two types of rack materials available:


Stainless steel:

- > Stable and robust for longevity and higher capacity
- > Fit freezer boxes with footprint up to 136 x 136 mm (5.4 x 5.4 in) for high flexibility
- > Optimized cooling conditions by special design for improved sample safety

Aluminum:


- > For freezer boxes up to 133 x 133 mm (5.2 x 5.2 in) footprint
- > Non-corrosive anodized aluminum for long-term use

- A. Stainless steel drawer rack
- B. Stainless steel side access rack
- C. Stainless steel rack for chest freezer
- D. Stainless steel rack for chest freezer with locking rod

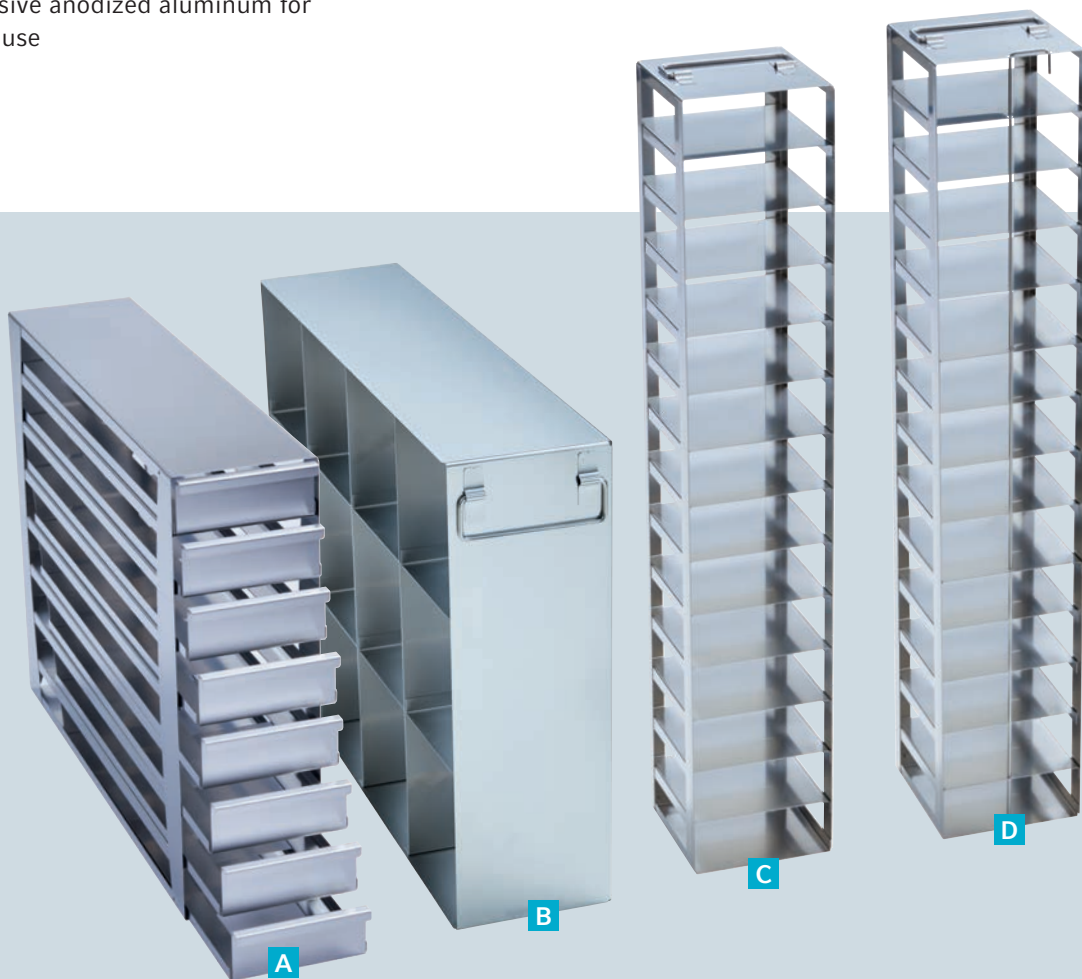

> Please check the rack brochure for easy rack selection or use the rack configuration tool at: <https://econfig.bio-itech.de>
 > Explore [The Eppendorf ULT Freezer Racking System Brochure](#) for more information



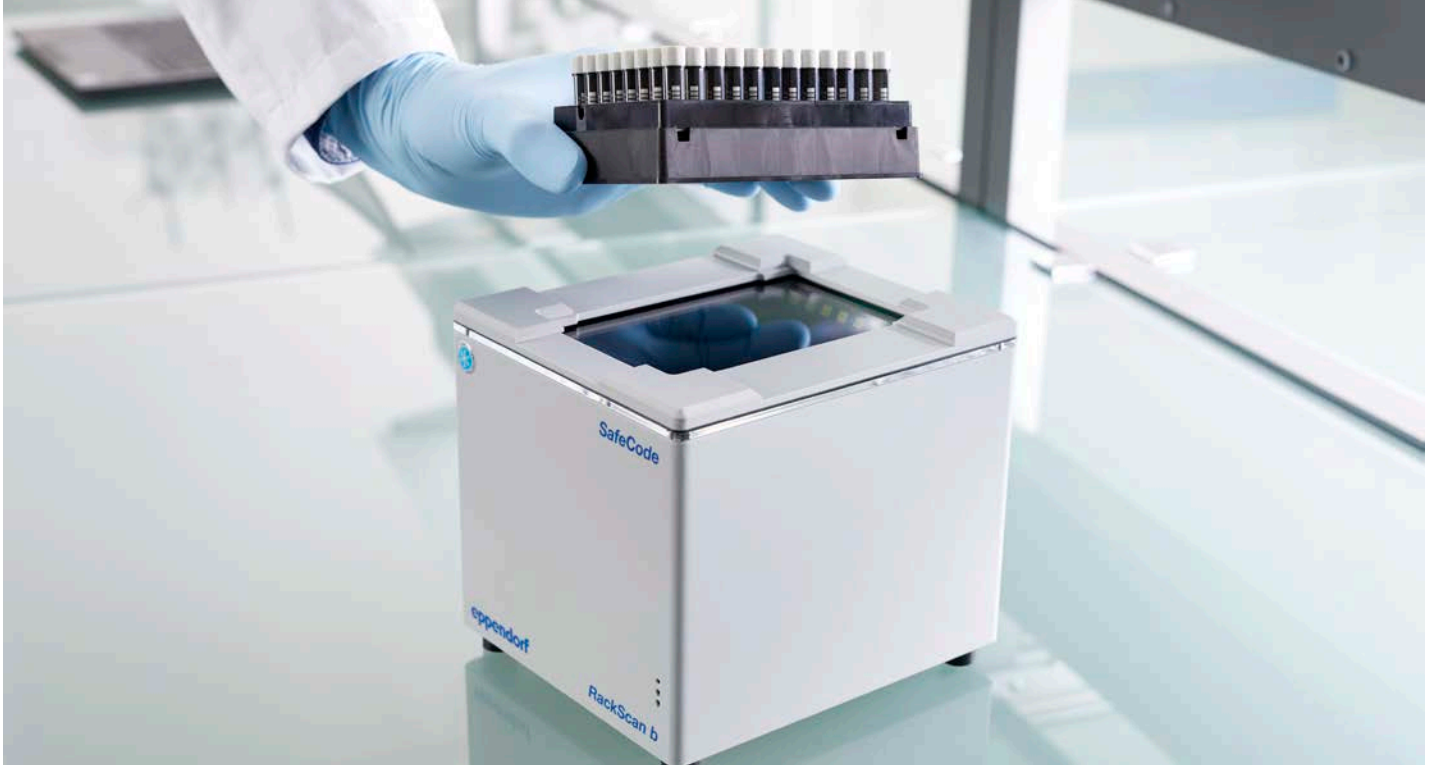

Safe Sample Identification Needed?

Secure your samples within the Eppendorf SafeCode System

CryoStorage Vials

The SafeCode system for Eppendorf CryoStorage Vials is based on 3-level coding to enable safe sample identification. Up to 96 pre-labeled vials are pre-racked and pre-capped for convenient usage.

- > Pre-labeled barcoded vessels with up to a 3-level-coding in different formats enable you to improve your processes
- > The ECC 200 error correction function enables safe code reading with up to 30% damaged code
- > Tube volumes from 0.5 mL – 4.0 mL provide a broad range of storage applications, available off-the-shelf for direct availability
- > Manufactured from high-quality polypropylene in a class VIII clean room for high purity level (Eppendorf Quality® or sterile purity)
- > Manage your barcoded samples with sample management software like eLabInventory
- > Receive further documentation about the vial from the Eppendorf DataPort: automatically downloaded by eLabInventory or via manual download

Scanning more than one in parallel?

The Eppendorf RackScan provides you different solutions to scan samples tubes with datamatrix (2D codes). The Eppendorf RackScan b scans data from the bottom of vessels. The instrument can be equipped with the „s“ accessory to read side-labeled vessel codes.

Data can be exported in different formats to downstream software applications like the eLabInventory software.




> www.eppendorf.com/SafeCode



Conical tubes

The tube range covers different volumes from 5 mL to 50.0 mL for a broad range of applications at the bench. This set of vessels expands your safe sample identification solution from the existing SafeCode CryoStorage Vials at -80 °C storage towards daily work at the bench.

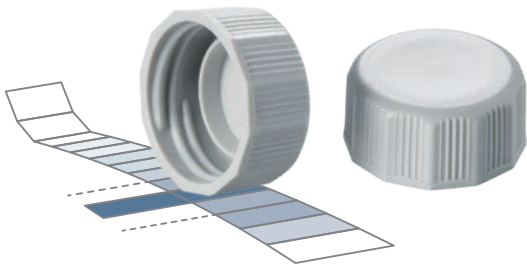

Cap with flattened and grooved sides

- > Supports user-friendly ergonomic handling
- > Enables stable and upright positioning
- > Minimizes contamination risk



Eppendorf DataPort



Do you require detailed documentation for lab certifications? Crucial data about your tubes, like **Lot Number**, **Order Number**, **Certificates**, etc., are available via the **Eppendorf DataPort website**. Export manually into databases or automatically to eLabInventory or eLabJournal®.

[Learn More!](#)



A system to meet your needs

The SafeCode System is designed to meet your needs at your convenience. Save time and effort with the correct tubes, apps, software, and storage solutions to make sure your valuable work ends up getting into publications and manufacturing, not the waste bin.

Lost in Samples? Stay Organized!

Over many years, scientists collect hundreds and thousands of samples – samples that are the results of years of hard work – samples of high value. When storing these, it is vital to keep them safe and accessible. Accurate sample labeling and proper storage in freezer boxes is only the beginning, you still need to maintain accurate records of the treasures in your freezer. Many scientists still use spreadsheets or even paper-based lists to keep track of their storage. While these solutions may work most of the time, they are susceptible to human error, damage, and misplacement.



Make science easier, make your life easier

Following the PhysioCare Concept, we recommend to use an intuitive, reliable sample management tool. Take advantage of modern solutions to make the lab experience more comfortable and more efficient.


eLabInventory

Keep track of your samples with eLabInventory


eLabInventory is a sample management software provided by Bio-ITech, an Eppendorf company. The tool organizes any item in the laboratory inventory, including specimens, materials, samples, and chemicals. Store samples in self-configured storage units such as freezers or refrigerators. eLabInventory is flexible and fully configurable to work in your laboratory.

- > Fully configurable to fit any type of laboratory
- > Stores any type of samples, specimens, and materials for high flexibility
- > Works on any device, delivered via the cloud or on-site
- > Intuitive user interface includes visual inventory browsing for easy handling
- > Barcode labeling and scanning for high-throughput
- > Import/export to other formats for data flexibility
- > Track and tracing for audit trail
- > GLP compliant sample management software
- > 21 CFR part 11 compliance for documentation security

Start a free 30-day trial!
Visit www.elabnext.com/eppendorf

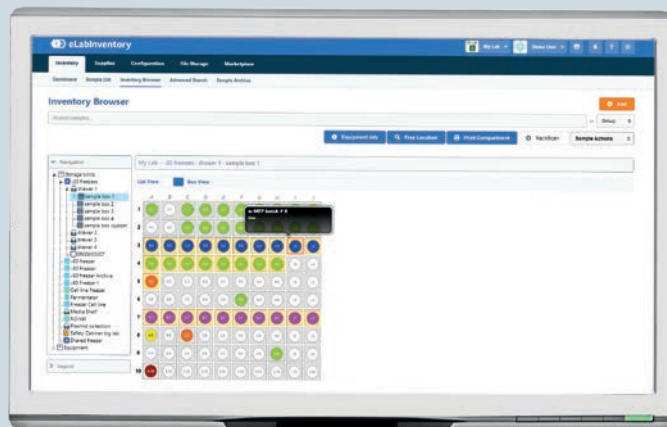

> For more information about ordering, visit:
www.elabnext.com/eppendorf

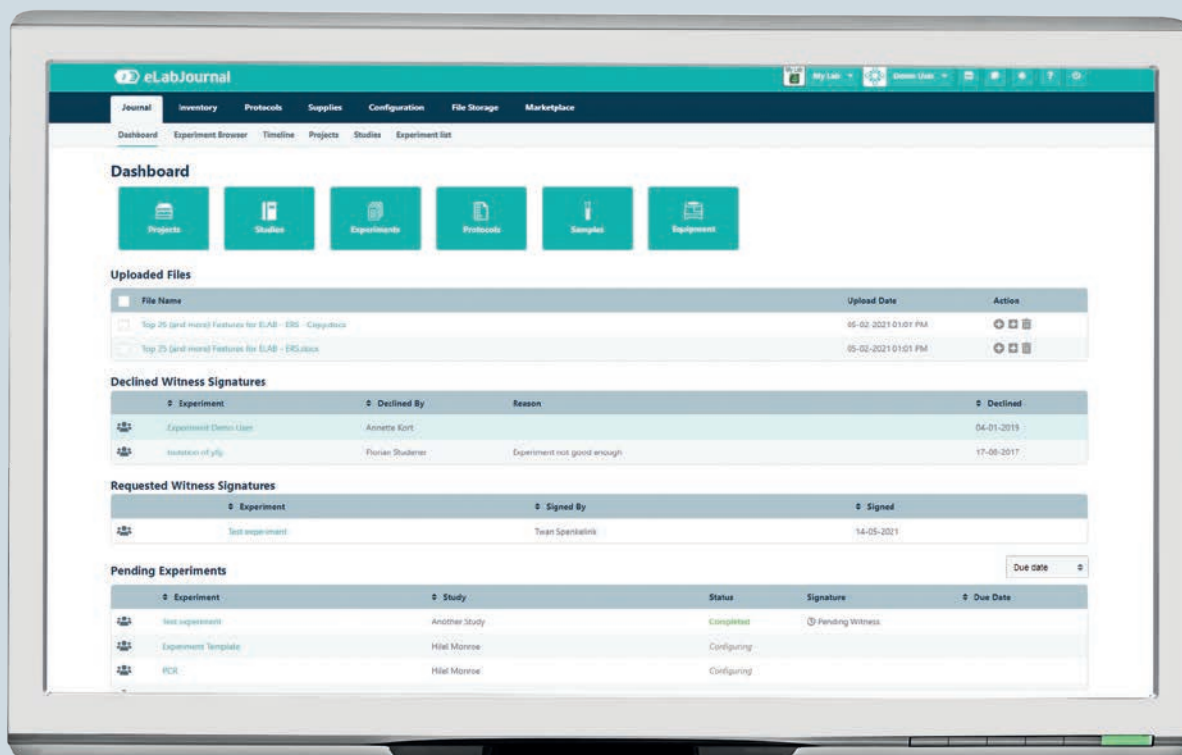

Start a free 30-day trial!
 Visit www.labnext.com/eppendorf

eLabJournal®

Improve your efficiency when documenting research and tracking research data with eLabJournal

The eLabJournal Electronic Lab Notebook offers an intuitive and flexible solution to document research and track research data. It improves efficiency when documenting, organizing, searching, and archiving collected data. With the free eLabJournal add-ons, you can extend the functionality of eLabJournal to a fully integrated Laboratory Information Management System. Designed by scientists for scientists.

- > Easily set up and apply workflows for lab routines using experiment templates
- > Use protocols as a template for your experiments and easily adjust the protocol according to your experiment settings using dynamic fields
- > Upload and link any type of file and link it to your experiment for safe documentation
- > eLabJournal has various add-ons to provide seamless integration and add great value to your account
- > To keep your data safe, we host our products in three data centers in different geographic locations that replicate and load-balance data in real time

> For more information, visit
www.elabnext.com/eppendorf



Your New ULT Has a History

Our future has history: Despite the more than 35 years of ULT experience with New Brunswick™, every new freezer project at Eppendorf begins with what we learn from before. Let's have a look at the CryoCube F740 series: The goal for R&D was to take what was great about our Innova® U725/-G ULT freezers and make it even better.



Yours
> Instrument is set-up in your lab




Delivery
> Market launch



Production
> Production of zero lot



Prototyping and beta-testing
> Prototyping
> Milestones for R&D
> Beta-tests with customers in the field
> Feedback from testing to be integrated in final R&D output

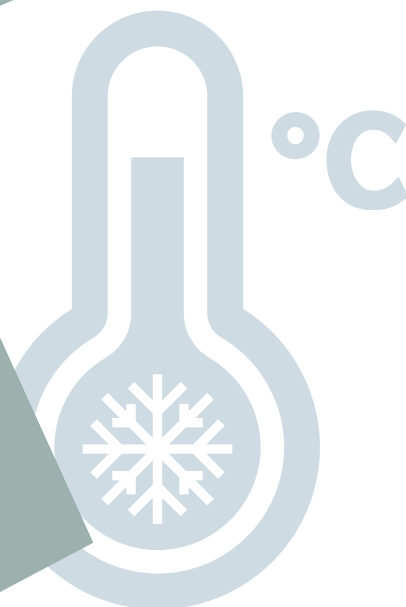


CAD
> Detailed specifications
> Precision engineering



Discussions and Learnings
> Internal findings: Interactive R&D with sales and production
> External findings: Feedback about existing freezers from users in the field, users of single units, as well as global pharma accounts and biobanks



Eppendorf CryoCube® F740hi

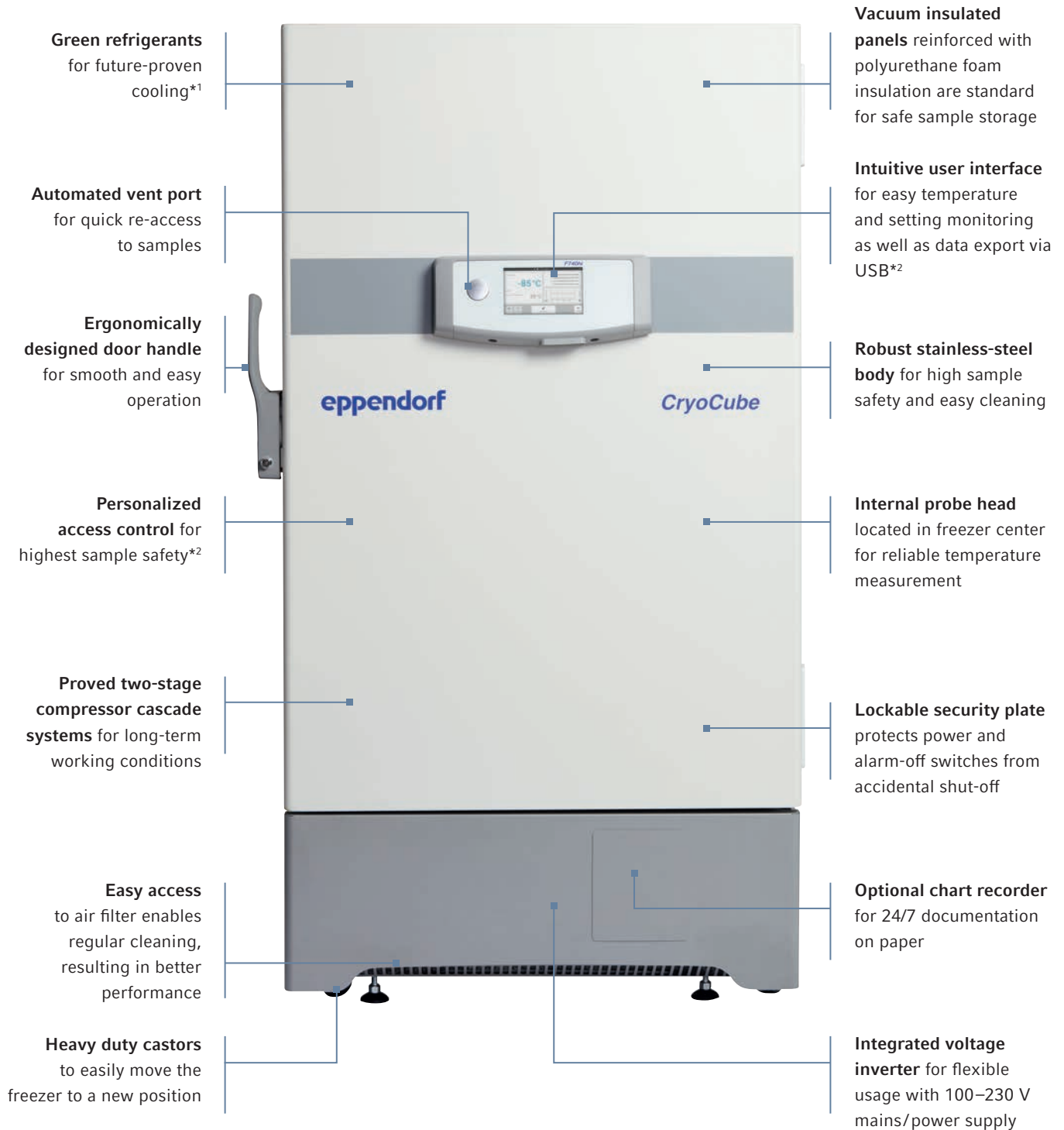

*¹ Available on "h" series models.

*² Available on "i" series models.



> **Ergonomic entry:** Comfortable door handle for easy opening with a mechanical lock



> **Ergonomic entry:** Comfortable door handle for easy opening
 > **Controlled access:** Personalized electronic access codes*², when necessary


Classic Dimensions, Modern Details




> **Get the maximum in:** 3-compartment set-ups with MAX racks use every empty spot*¹



> **Easy access:** Quickly open the bottom-mounted air filter for easy cleaning and maintenance*¹
 > **Simple solution:** No tools required to access the filter

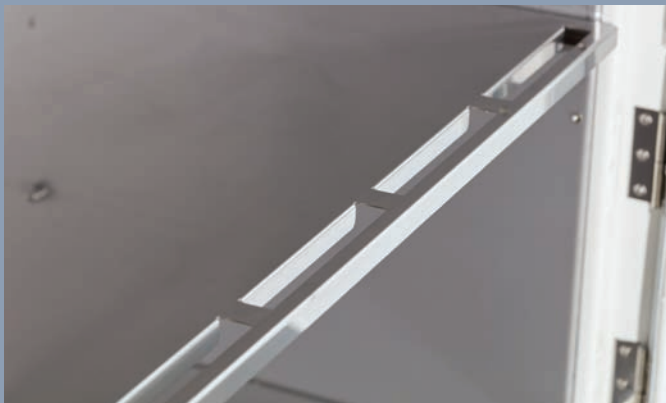

> **Enhanced shelf design:** Integrated airways for fast and consistent air circulation*¹
 > **Keep your fingers safe:** Rounded metal rims at shelves

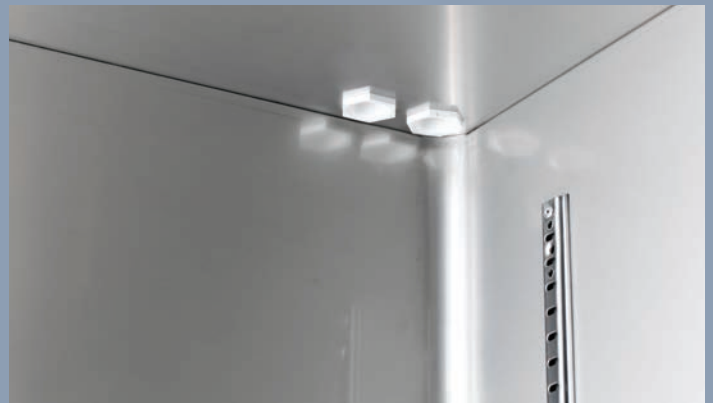

> **External probes or back-ups:** Two access ports in the upper corner
 > **Easy cleaning:** Stainless steel inner chamber

*¹ Available on CryoCube F740 series models.
 *² Available on CryoCube F740 "i" series models.
 *³ Available on CryoCube series models.



> **Concentrate on the essentials:** Easy settings



> **If you are in a hurry:** Automatic vent port enables fast sample access^{*3}



> **Decide what you need:** 3-compartment/shelves or 5-compartment/shelves with 3 inner doors^{*1}



> **Twins except the door handle^{*1}:** Decide what fits best in your lab, a left-handed door handle or a right-handed door handle



> **Make life easier:** Inner doors can easily be removed for full access to compartments during cleaning of freezer



> **If you prefer privacy:** The CryoCube F101h – a real ULT freezer in one small package below the bench

^{*1} Available on CryoCube F740 series models.



Model	CryoCube® F740hi, 3 c./ CryoCube® F740hi, 5 c.	CryoCube® F740hiw, 3 c./ CryoCube® F740hiw, 5 c.	CryoCube® F740h, 3 c./ CryoCube® F740h, 5 c.
Insulation	Vacuum insulation paneling / polyurethane foam	Vacuum insulation paneling / polyurethane foam	Vacuum insulation paneling / polyurethane foam
Wall Thickness	Type B	Type B	Type B
Capacity	740 L / 26.1 ft ³	740 L / 26.1 ft ³	740 L / 26.1 ft ³
No. of internal doors	3	3	3
No. of compartments	3/5	3/5	3/5
Max. racks per freezer	18/30	18/30	18/30
Box capacity per rack			
53 mm (2 in) tall boxes	32/16	32/16	32/16
76 mm (3 in) tall boxes	20/8	20/8	20/8
102 mm (4 in) tall boxes	12	12	12
Box capacity per freezer			
53 mm (2 in) tall boxes	576/528	576/528	576/528
76 mm (3 in) tall boxes	360/288	360/288	360/288
102 mm (4 in) tall boxes	288/264	288/264	288/264
Sample capacity			
53 mm (2 in)	57,600/52,800	57,600/52,800	57,600/52,800
Dimensions (W × D × H) ¹⁾²⁾³⁾	111 × 98 × 197.3 cm	111 × 98 × 197.3 cm	111 × 91.5 × 197.3 cm
Dimensions internal (W × D × H)	86.5 × 62.1 × 139 cm	86.5 × 62.1 × 139 cm	86.5 × 62.1 × 139 cm
Weight w/o accessories	315 kg/333 kg	320 kg/328 kg	308 kg/317 kg
Shipping weight	364 kg/382 kg	369 kg/377 kg	357 kg/ 366 kg
Pull-down time to -80 °C (230 V) ⁴⁾	4 h 10 min	4 h 10 min	4 h 10 min
Noise level	41.3 dB	41.3 dB	47.8 dB
Heat output (230 V)	438 W	442 W	492 W
Power supply	100-230 V, 50/60 Hz	100-230 V, 50/60 Hz	100-230 V, 50/60 Hz
Power consumption	10.5 kWh/day	10.6 kWh/h	11.8 kWh/day
Eppendorf freezer (230/115 V) ⁵⁾	11.6 kWh/day	11.7 kWh/day	13.0 kWh/day
Interior	Stainless steel	Stainless steel	Stainless steel
Password protection of setpoints	yes	yes	yes
Extras	Automatic vent port	Automatic vent port	Automatic vent port
Cooling liquid	R290 / R170	R290 / R170	R290 / R170
Connectivity	VisioNize touch enabled	VisioNize touch enabled	VisioNize box

¹⁾ Optional CO₂/LN₂ backup systems add 8.65 cm/3.5 in to height.

²⁾ For usage of handles and hinges, add 80 mm to width of upright freezers and 110 mm to the depth of chest freezers.

³⁾ Door open adds up to 15 cm.

⁴⁾ Empty freezer with shelves fitted, upright freezers only, pull down from 20 – 22 °C ambient conditions.

⁵⁾ Empty freezer with shelves fitted, upright freezers only, set point -80 °C, 20 – 22 °C ambient conditions.



> For more information go to www.eppendorf.com/freezers

**CryoCube® F570h**Vacuum insulation paneling /
polyurethane foam

Type C

570 L / 20.1 ft³

5

5

25

16

8

8

400

200

200

40,000

110.3 x 89.8 x 195 cm (43.4 x 35.4 x 76.8 in)

76.5 x 57.5 x 126.5 cm (30.1 x 22.6 x 49.8 in)

341 kg (752 lbs)

388 kg (855 lbs)

3 h 30 min

57.7 dB

308 W

230 V (50 Hz), 115 V (60 Hz)

7.4 KWh/day

7.1 KWh/day

Stainless steel

yes

Automatic vent port

R290 / R170

VisioNize box

CryoCube® F570hwVacuum insulation paneling /
polyurethane foam

Type C

570 L / 20.1 ft³

5

5

25

16

8

8

400

200

200

40,000

110.3 x 89.8 x 195 cm (43.4 x 35.4 x 76.8 in)

76.5 x 57.5 x 126.5 cm (30.1 x 22.6 x 49.8 in)

341 kg (752 lbs)

388 kg (855 lbs)

4 h 20 min

56.0 dBA

325 W

230 V (50 Hz)

8.0 KWh/day

-

Stainless steel

yes

Automatic vent port, water-cooled

R290 / R170

VisioNize box

CryoCube® F570n

Polyurethane foam

Type A

570 L / 20.1 ft³

5

5

25

16

8

8

400

200

200

40,000

110.3 x 89.8 x 195 cm (43.4 x 35.4 x 76.8 in)

76.5 x 57.5 x 126.5 cm (30.1 x 22.6 x 49.8 in)

281 kg (619 lbs)

328 kg (723 lbs)

3 h 50 min

57.5 dB

346 W

230 V (50 Hz), 208–230 V (60 Hz), 115 V (60 Hz)

8.3 KWh/ day (208 V: 8.5 KWh/ day)

8.8 KWh/ day

Stainless steel

yes

Automatic vent port

R290 / R170

VisioNize box





Model	Innova® U535	CryoCube® F440h
Insulation	Vacuum insulation paneling / polyurethane foam	Vacuum insulation paneling / polyurethane foam
Wall Thickness	Type B	Type C
Capacity	535 L / 18.9 ft ³	440 L / 15.5 ft ³
No. of internal doors	3	5
No. of compartments	3	5
Max. racks per freezer	12	20
Box capacity per rack		
53 mm (2 in) tall boxes	28	16
76 mm (3 in) tall boxes	20	8
102 mm (4 in) tall boxes	12	8
Box capacity per freezer		
53 mm (2 in) tall boxes	336	320
76 mm (3 in) tall boxes	240	160
102 mm (4 in) tall boxes	144	160
Sample capacity		
53 mm (2 in)	33,600	32,000
Dimensions (W × D × H) ¹⁾²⁾³⁾	80 × 86.7 × 195 cm (31.5 × 34.1 × 76.8 in)	95.5 × 89.8 × 195.0 cm (37.6 × 35.4 × 76.8 in)
Dimensions internal (W × D × H)	64 × 61.5 × 136.5 cm (25.2 × 24.2 × 53.7 in)	61.6 × 57.5 × 126.5 cm (24.3 × 35.4 × 49.8 in)
Weight w/o accessories	250 kg (551 lbs)	286 kg (631 lbs)
Shipping weight	296 kg (653 lbs)	327 kg (721 lbs)
Pull-down time to -80 °C (230 V) ⁴⁾	4 h	3 h 20 min
Noise level	56 dBA	51.5 dBA
Heat output (230 V)	550 W	283 W
Power supply	230 V (50 Hz), 208-230 V (60 Hz), 115 V (60 Hz)	230 V (50 Hz), 115 V (60 Hz)
Power consumption	13.2 kWh/day	6.8 kWh/day
Eppendorf freezer (230/115 V) ⁵⁾	13.2 kWh/day	6.5 kWh/day
Interior	Stainless steel	Stainless steel
Password protection of setpoints	yes	yes
Extras	-	-
Cooling liquid	R404A / R508B	R290 / R170
Connectivity	VisioNize box	VisioNize box

¹⁾ Optional CO₂/LN₂ backup systems add 8.65 cm/3.5 in to height.

²⁾ To allow for handles and hinges, add 80 mm to width of upright freezers and 110 mm to the depth of chest freezers.

³⁾ Door open adds up to 15 cm.

⁴⁾ Empty freezer with shelves fitted, upright freezers only, pull down from 20 – 22 °C ambient conditions.

⁵⁾ Empty freezer with shelves fitted, upright freezers only, set point -80 °C, 20 – 22 °C ambient conditions.












> For more information go to www.eppendorf.com/freezers



CryoCube® F440n	CryoCube® FC660h	CryoCube® F101h
Polyurethane foam	Vacuum insulation paneling / polyurethane foam	Vacuum insulation paneling / polyurethane foam
Type A	Type C	Type B
440 L / 15.5 ft ³	660 L / 23.3 ft ³	101 L / 3.6 ft ³
5		2
5	N/A	2
20	40	6
16	13	10
8	9	6
8	5	4
320	520	60
160	360	36
160	280	24
32,000	52,000	6,000
95.5 x 89.8 x 195.0 cm (37.6 x 35.4 x 76.8 in)	207 x 94.5 x 107.5 cm (81.5 x 37.2 x 42.3 in)	93.6 x 62.9 x 83.4 cm (36.9 x 24.8 x 32.8 in)
61.6 x 57.5 x 126.5 cm (24.3 x 35.4 x 49.8 in)	147 x 59 x 76 cm (57.9 x 23.2 x 29.9 in)	48 x 33 x 64 cm (18.9 x 13 x 25.2 in)
267 kg (589 lbs)	304 kg (670 lbs)	123 kg (271 lbs)
308 kg (679 lbs)	379 kg (836 lbs)	153 kg (337 lbs)
3 h 35 min	2 h 50 min	2 h
51.5 dB	58.3 dBA	39.7/43.2 dBA
346 W	338 W	196 W
230 V (50 Hz), 208-230 V (60 Hz), 115 V (60 Hz)	230 V (50 Hz), 208-230 V (60 Hz)	230 V (50 Hz), 115 V (60 Hz)
8.3 KWh/day	8.1 KWh/day	4.7 KWh/day
8.5 KWh/day	9.5 KWh/day	5.1 KWh/day
Stainless steel	Stainless steel	Stainless steel
yes	yes	yes
-	30 mm thick inner lids for temp uniformity	Fits below bench
R290 / R170	R290 / R170 (R404A / R508B in USA)	R290 / R170
VisioNize box	VisioNize box	VisioNize box







Freezer Ordering information

Model	Right or Left Door Handle	Inner Compartments	Voltage*	Order No.
CryoCube® F740h	 Left	5	115 – 208 V	F740540035**
			230 V	F740340031
		3	115 – 208 V	F740540015**
	 Right	5	115 – 208 V	F740540045**
			230 V	F740340041
		3	115 – 208 V	F740540025**
			230 V	F740340021
CryoCube® F740hi	 Left	5	115 – 208 V	F740520035**
			230 V	F740320031
		3	115 – 208 V	F740520015**
	 Right	5	115 – 208 V	F740520045**
			230 V	F740320041
		3	115 – 208 V	F740520025**
			230 V	F740320021
CryoCube® F740hiw Water-cooled	 Left	5	115 – 208 V	F740520135**
			230 V	F740320131
		3	115 – 208 V	F740520115**
			230 V	F740320111
CryoCube® F570n	 Left	5	115 V	F571200035
			208 V	F571400035
			230 V	F571300031
CryoCube® F570h	 Left	5	115 V	F571240035
			230 V	F571340031
CryoCube® F570hw Water-cooled	 Left	5	230 V	F571340131
Innova® U535	 Left		115 V	U9430-0000
			208 V	U9430-0002
			230 V	U9430-0001

*230 V are 50–60 Hz with European (Schuko) plug, further plug versions available
 **CryoCube F740 series (115/208 V): Equipped with both cables by standard

Freezer Ordering information

Model	Right or Left Door Handle	Inner Compartments	Voltage*	Order No.
CryoCube® F440h 	Left	5	115 V	F440240035
			230 V	F440340031
CryoCube® F440n 	Left	5	115 V	F440200035
			208 V	F440400035
			230 V	F440300031
CryoCube® F101h 	Right	2	115 V	F101240065
			230 V	F101340061
CryoCube® FC660h 	-	-	208 V	F660420005
			230 V	F660320001

*230 V are 50-60 Hz with European (Schuko) plug, further plug versions available

VisioNize® Ordering Information







Description	Order no.
VisioNize® box, with mains/power adapter	1006 000 009
VisioNize® box adapter for freezers, CryoCube® and Innova® models	1006 075 009
VisioNize® box cable, Ethernet, length 5 m	1006 073 006
RS-485 interface kit, for Premium/Innova®/HEF®/CryoCube® F570 series (except Innova® U360)	F652 999 007



Discover the entire SafeCode System, including CryoStorage Vials, Eppendorf RackScan, SafeCode Barcoded Tubes & Plates, eLabInventory, eLabJournal, & much more! www.eppendorf.com/SafeCode



Ordering Information

Description	Order no.	
Eppendorf Tubes® 5.0 mL , Sterile, SafeCode, 200 tubes	0030 113.500	
Eppendorf Conical Tubes 15 mL , Sterile, SafeCode, 500 tubes	0030 113.489	
Eppendorf Conical Tubes 50 mL , Sterile, SafeCode, 500 tubes	0030 113.497	
Eppendorf CryoStorage Vials , SafeCode portfolio	See online	
DWP 96/2000 , SafeCode, 20x (5 bags x 4 plates)	0030 113.527	
DWP 96/1000 , SafeCode, 20x (5 bags x 4 plates)	0030 113.535	
DWP 96/500 , SafeCode, 40x (5 bags x 8 plates)	0030 113.543	
DWP 384/200 , SafeCode, 40x (5 bags x 8 plates)	0030 113.551	
PCR twin.tec 96 , skirted, SafeCode, 25x	0030 113.560	
PCR twin.tec 384 , skirted, SafeCode, 25x	0030 113.578	
MTP 96/V clear , PCR clean, SafeCode, 80x (5 bags x 16 plates)	0030 113.586	
Eppendorf RackScan b , scanner (bottom reading) for 2D-coded tubes in SBS racks and cryostorage boxes, incl. RackScan Analyzer software, 230 V/50 - 60 Hz	6002 000 006	
Eppendorf RackScan s , code scanner for side reading of SBS racks and cryostorage boxes, accessory for the RackScan b	6002 000.014	
	Eppendorf CryoStorage Vial 0.5 mL , pre-capped, pre-racked, 10x96 vials, 2D SafeCode	0030 079.400
	Eppendorf CryoStorage Vial 0.5 mL , pre-capped, pre-racked, 10x96 vials, 2D SafeCode, sterile	0030 079.450
	Eppendorf CryoStorage Vial 1.0 mL , pre-capped, pre-racked, 10x96 vials, 2D SafeCode	0030 079.418
	Eppendorf CryoStorage Vial 1.0 mL , pre-capped, pre-racked, 10x96 vials, 2D SafeCode, sterile	0030 079.469
	Eppendorf CryoStorage Vial 1.5 mL , pre-capped, pre-racked, 10x48 vials, 2D SafeCode	0030 079.426
	Eppendorf CryoStorage Vial 1.5 mL , pre-capped, pre-racked, 10x48 vials, 2D SafeCode, sterile	0030 079.477
	Eppendorf CryoStorage Vial 2.0 mL , pre-capped, pre-racked, 10x48 vials, 2D SafeCode	0030 079.434
	Eppendorf CryoStorage Vial 2.0 mL , pre-capped, pre-racked, 10x48 vials, 2D SafeCode, sterile	0030 079.485
	Eppendorf CryoStorage Vial 4.0 mL , pre-capped, pre-racked, 10x48 vials, 2D SafeCode	0030 079.442
	Eppendorf CryoStorage Vial 4.0 mL , pre-capped, pre-racked, 10x48 vials, 2D SafeCode, sterile	0030 079.493

Start a free 30-day trial!
Visit www.elabnext.com/eppendorf



eLabInventory

- > Stores any type of samples, specimens, and materials for high flexibility
- > Intuitive user interface includes visual inventory browsing for easy handling
- > Barcode labeling and scanning for high-throughput
- > Fully customizable data records for samples, storage locations, and equipment
- > Import/export to other formats for data flexibility
- > 21 CFR part 11 compliant for documentation security

Ordering Information

Description	Industrial Order no.	Academic Order no.
eLabInventory Cloud, 3-seat year license	BIT0140503	BIT0140501
eLabInventory Cloud, 5-seat year license	BIT0140103	BIT0140101
eLabInventory Cloud, 10-seat year license	BIT0140203	BIT0140201
eLabInventory Cloud, 25-seat year license	BIT0140303	BIT0140301
eLabInventory Cloud, 50-seat year license	BIT0140403	BIT0140401

System maintenance, software updates, and support are included in the offered licenses.

Start a free 30-day trial!
Visit www.elabnext.com/eppendorf



eLabJournal

- > Intuitive user interface for easy handling
- > Use protocols as a template for your experiments and easily adjust the protocol according to your experiment setting using dynamic fields to save time
- > Upload and link any type of file and link it to your experiment for safe documentation
- > eLabJournal has various add-ons to provide seamless integration and expand your possibilities
- > 21 CFR part 11 compliant for documentation security

Ordering Information

Description	Industrial Order no.	Academic Order no.
eLabJournal® Cloud, 3-seat year license	BIT0100503	BIT0100501
eLabJournal® Cloud, 5-seat year license	BIT0100103	BIT0100101
eLabJournal® Cloud, 10-seat year license	BIT0100203	BIT0100201
eLabJournal® Cloud, 25-seat year license	BIT0100303	BIT0100301
eLabJournal® Cloud, 50-seat year license	BIT0100403	BIT0100401

System maintenance, software updates, and support are included in the offered licenses.



CO₂ and LN₂ Back-up System

- > Equipped with a battery back-up to temporarily protect the contents of your freezer in case of power failure
- > Available for either liquid CO₂ and LN₂. Liquid CO₂ can maintain temperatures between -50 °C and -70 °C, where as LN₂ can go down to -85 °C
- > Both systems can be factory installed with your freezer order or retrofitted on location
- > Info: CO₂ and LN₂ are not interchangeable systems

Ordering Information

Description	Order no.
Innova® U101; U360; U535; U725/ -G; C585; C760; CryoCube® F101h; F740 series	
CO ₂ back-up system, 100 V/50 – 60 Hz, Innova® (»narrow« version)	F652999005
CO ₂ back-up system, 120 – 220 V/60 Hz, Innova® (»narrow« version)	U9043-0002
CO ₂ back-up system, 230 V/50 Hz, Innova® (»narrow« version)	U9043-0004
LN ₂ back-up system, 100 V/50 – 60 Hz, Innova® (»narrow« version)	F652999006
LN ₂ back-up system, 120 – 220 V/60 Hz, Innova® (»narrow« version)	U9044-0002
LN ₂ back-up system, 230 V/50 Hz, Innova® (»narrow« version)	U9044-0004
CryoCube® F440 series, F570 series, FC660 series; Premium U410; HEF U410	
CO ₂ back-up system, 120 – 220 V/60 Hz, (»broad« version)	U9043-0006
CO ₂ back-up system, 230 V/50 Hz, (»broad« version)	U9043-0008
LN ₂ back-up system, 120 – 220 V/60 Hz, (»broad« version)	U9044-0006
LN ₂ back-up system, 230 V/50 Hz, (»broad« version)	U9044-0008

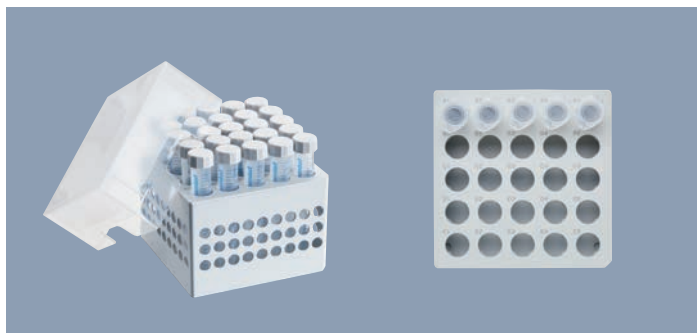


Chart Recorder

- > Includes chart recorder, RTD probe, keys, chart recorder paper, and chart recorder pens
- > Front mounted to your freezer for easy access and can be fitted to all freezer models
- > User configurable built in high and low temperature alarms
- > Can be factory installed with your freezer order or retrofitted on location
- > Selectable temperature range (0 to -50 °C and -50 to -100 °C) and rotation speed
- > Type I: Battery powered, will provide a minimum of eighteen months of operation

Ordering Information

Description	Order no.
Chart Recorder (Type I)	P0625-2100
Chart Recorder pens (Type I), 3 pk	K0660-0051
Chart Recorder paper (Type I), 0 to -50 °C	P0625-2111
Chart Recorder paper (Type I), -50 to -100 °C	P0625-2110



Freezer Storage Boxes:

- > Outer footprint dimensions (133 mm) are compatible with common freezer rack systems
- > Inner grid variants for all typical tube formats: cryogenic tubes, microcentrifuge tubes, conical tubes 15 and 50 mL, and other laboratory vessels

Ordering information

Description	Order no.
Storage Box 10 × 10 , for 100 cryogenic tubes w. int. thread, 3 pcs., height 52.8 mm, 2 inch, polypropylene, for freezing to -86 °C, autoclavable, with lid and alphanumeric code	0030 140.508
Storage Box 9 × 9 , for 81 screw cap (cryog.) tubes 1-2 mL, 3 pcs., height 52.8 mm, 2 inch, polypropylene, for freezing to -86 °C, autoclavable, with lid and alphanumeric code	0030 140.516
Storage Box 9 × 9 , for 81 screw cap (cryog.) tubes 3 mL, 2 pcs., height 76.2 mm, 3 inch, polypropylene, for freezing to -86 °C, autoclavable, with lid and alphanumeric code	0030 140.540
Storage Box 9 × 9 , for 81 screw cap (cryog.) tubes 4-5 mL, 2 pcs., height 101.6 mm, 4 inch, polypropylene, for freezing to -86 °C, autoclavable, with lid and alphanumeric code	0030 140.567
Storage Box 8 × 8 , for 64 tubes 1-2 mL, 3 pcs., height 52.8 mm, 2 inch, polypropylene, for freezing to -86 °C, autoclavable, with lid and alphanumeric code	0030 140.524
Storage Box 5 × 5 , for 25 tubes 5 mL, 4 pcs., height 63.5 mm, 2.5 inch, polypropylene, for freezing to -86 °C, autoclavable, with lid and alphanumeric code	0030 140.532
Storage Box 5 × 5 , for 25 tubes 5 mL screw cap, 2 pcs., height 76.2 mm, 3 inch, polypropylene, for freezing to -86 °C, autoclavable, with lid and alphanumeric code	0030 140.613
Storage Box 5 × 5 , for 25 tubes 15 mL, 2 pcs., height 127 mm, 5 inch, polypropylene, for freezing to -86 °C, autoclavable, with lid and alphanumeric code	0030 140.583
Storage Box 3 × 3 , for 9 tubes 25 mL, 2 pcs., height 89 mm, 3.5 inch, polypropylene, for freezing to -86 °C, autoclavable, with lid and alphanumeric code	0030 140.729
Storage Box 3 × 3 , for 9 tubes 50 mL and 4 tubes 15 mL, 2 pcs., height 127 mm, 5 inch, polypropylene, for freezing to -86 °C, autoclavable, with lid and alphanumeric code	0030 140.591

> Ordering information regarding racks is available online, via the online configurator, or within a separate brochure

Fitness for Your ULT Freezer



Take good care of your freezer

Taking care of your freezer can help ensure optimal storage of your samples for years to come.



Therefore, we strongly recommend that you perform the following actions on a regular basis:

Freezer maintenance		As needed	3 months	6 months	12 months
Cleaning door gaskets	Dry cloth, brush	■			
Remove snow/ thin ice	Cloth, brush	■			
Vacuum relief port	Cloth	■			
Filter	Vacuum	■	■		
Condenser cleaning	Vacuum			■	
Defrost and cleaning	Power off, dry				■

Ordering information

Description	Order no.
Freezer air filter	
Innova® U101	K0200-0506
Innova® U360	P0625-1110
CryoCube® F101h	K0200-0506
CryoCube® F740 series	F740850029
Non-listed upright models	K0200-0511
Chest Models	K0200-0516
Freezer extra shelf , includes four mounting clips	
Innova® U101	P0625-0180
Innova® U360	P0625-1180
Innova® U535	K0280-1034
Innova® U725 series	K0280-1036
U410 series	K0280-1031
CryoCube® F101h (2023 Update)	F652999009
CryoCube® F570 series	K0280-1030
CryoCube® F440 series	F651999580
CryoCube® F740 series	F740850125
Stainless Steel Shelf Clips, pack of 4	K0280-0550
Inner lids for Chest Freezers	
Innova® C585	K0160-0777
Innova® C760	K0160-0776
CryoCube® FC660 series	K0280-1043

Ordering information

Description	Order no.
Additional	
Battery, 6V, 2.8 Ah, all models except CryoCube® F740 series, Innova® U360	K0380-0170
Battery (Li) for CryoCube® F740 series	F740850003
Battery for Innova® U360	P0625-0970
Alarm connector plug, for bms relay	K0380-0451
External Voltage Stabilizer, 120V/60 Hz	M1322-0000
External Voltage Stabilizer, 208-230V/60 Hz	M1322-0002
Innova® U101 stacking kit	P0625-0020
Padlock adapter kit, for upright models except CryoCube® models	P0625-1170
Padlock adapter for CryoCube F101h (2023 Update)	F652999010

Pit Stop Needed?

In addition to our extensive standard warranty, Eppendorf offers ULT Freezer Performance Plans to suit your lab's needs while enhancing your peace of mind. You can choose between essential function checks, preventive maintenance plans, or Installation and Operational (IQ/OQ) certification services. Protect your priceless samples!

Preventive Maintenance

The Eppendorf Ultra-low Temperature Freezer Performance Plans components include a choice of preventive maintenance protocols designed to check, validate and to promote continuous stable freezer performance, optimizing instrument conformity with programmed storage temperatures.

- > The **Device Check** includes a check of all fundamental functions of the product.
- > The **Preventive Maintenance** includes all preventive maintenance services in order to keep the instrument in accordance with the manufacturer specifications.
- > Our preventive maintenance offerings ensure that stored samples are protected through continuous refrigeration efficiency, reliability, and rapid recovery of storage temperatures to pre-set levels following door openings.

Certification Services

Installation Qualification (IQ) and Operational Qualification (OQ) certification services support your Quality Management requirements by providing you with qualified assurance that your ultra-low temperature freezer is functioning correctly, in accordance with manufacturer specifications.



eppendorf

Performance tested on

Model: _____

Serial no.: _____

Serviced by: _____

Service no.: _____

Next service: _____

Date: _____

epServices
for premium performance

»Eppendorf ULT Freezers: Your Samples are Safe with Us.«

View even more freezers, accessories, and options

Eppendorf ULT Freezers have established a reputation for outstanding sample quality preservation while minimizing energy consumption. Combined with rapid delivery and service, Eppendorf freezers continue to raise the bar.


> Explore all the cold storage solutions Eppendorf has to offer on our extensive web page and eShop



> Learn more about sample security:

www.eppendorf.com/freezers

> Have a detailed look at our entire ULT Freezer Portfolio using this QR Code

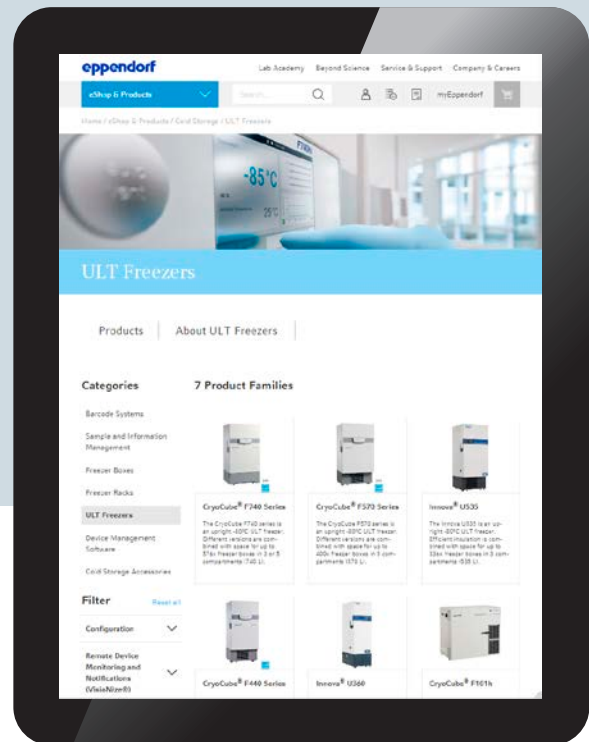

Your local distributor: www.eppendorf.com/contact
Eppendorf SE · Barkhausenweg 1 · 22339 Hamburg · Germany
eppendorf@eppendorf.com · www.eppendorf.com

www.eppendorf.com/freezers